

SafePass Fluid Pass Thru System

Sentinel's SafePass Fluid Pass Thru System provides users with an economical way to install tubing transfer lines between classified spaces. Easy to install and operate, SafePass provides a positive seal between suites that protects the integrity of the cleanroom environment. The engineered, precision machined Isolation Seal provides for aseptic transfer of fluid between rooms. SafePass's unique design prevents unwanted airflow between rooms providing complete isolation.

Sentinel's SafePass System is designed to accept all tubing brands and materials or custom tube sets to transfer fluid from suite to suite.


- Multi-directional installation based on process conditions and requirements
- No proprietary tubing required, use any tubing material and brand you desire
- ISO 7 customized tube sets available to your specifications, Gamma Irradiation available
- Available in multiport or single port configuration
- Prevents unwanted airflow between suites
- Pressure rated system - seals to 1psi
- Simple replacement of critical seals
- Designed for tubing sizes from 1/4" to 3/4" OD

Visit us at sentinelprocess.com or call 800.345.FLOW today!


Contact Sentinel Process Systems, Inc. to help with all your Process Component, Instrumentation, and Single-Use needs.



Fluid Pass Thru System

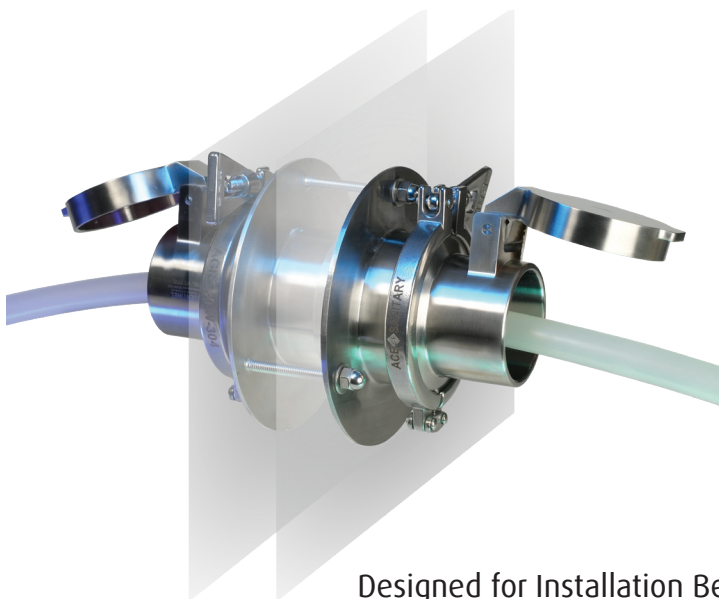

Designed for Installation Between Modular Wall Systems

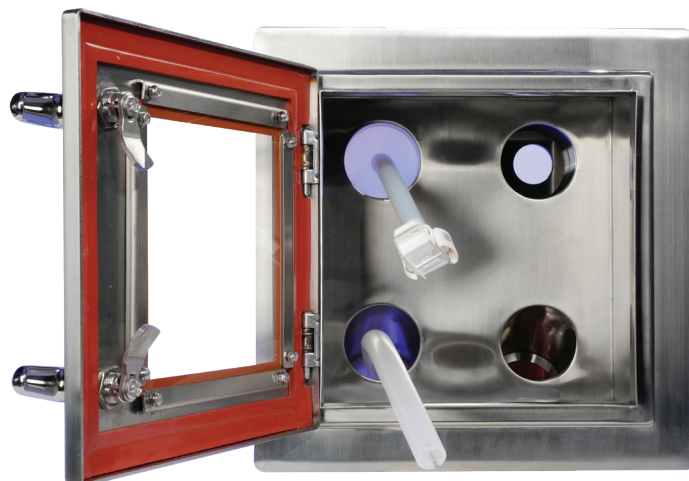

Multiport Cabinet Design Available with Sight Glass Port, Straight or Sweep Outlets



Delivering Innovative Process Technologies



800.345.FLOW
sentinelprocess.com