

# CUSTOMER APPROVAL

COMPANY:


APPROVED BY:

DATE:

NOTE: CUSTOMER SIGNATURE INDICATES APPROVAL OF DRAWING AND ACCEPTANCE OF ALL COMPLIANCE AND MATERIALS AS REPRESENTED ON THE DRAWING.

## PARTS LIST

ITEM	QTY	TITLE	PART NUMBER	DESCRIPTION
1	44.000 in	ADVANTAPURE PLATINUM CURED SILICONE	APST-0375-0625	3/8" I.D. x 5/8" O.D.
2	20.000 in	ADVANTAPURE HIGH PRESSURE PLATINUM CURED SILICONE	APHP-0375-0625	3/8" I.D. x 5/8" O.D.
3	1	FLWSENSOR	LFS06SUZSC1	LEVIFLOW ULTRASONIC SU FLOW METER, POLYPROPYLENE
4	1	QUATTROFLOW PUMP	QF15DISPP-3-EZ	MINI-CLAMP I/O, POLYPROPYLENE
5	1	SciLog Pressure Sensor II Series	206-252	3/8" HOSE BARB I/O
6	4	ASEPTIQUIK G CONNECTOR	AQG17006	3/8" HOSE BARB - BLUE PULL TAB
7	4	PLATINUM-CURED SILICONE GASKET		3/4" MINI-CLAMP
8	4	MINI TC CLAMP		1/2" to 3/4" TC, WHITE NYLON
9	4	TUBE FITTING WITH MINI T.C.		3/8 HOSE BARB X 3/4" T.C. - POLYPROPYLENE
10	2	TUBE FITTING - TEE		3/8" HOSE BARBS - POLYPROPYLENE
11	1	PINCH CLAMP		1/2"-5/8" O.D.
12	16	OETIKER CLAMP		5/8in OD 304SS

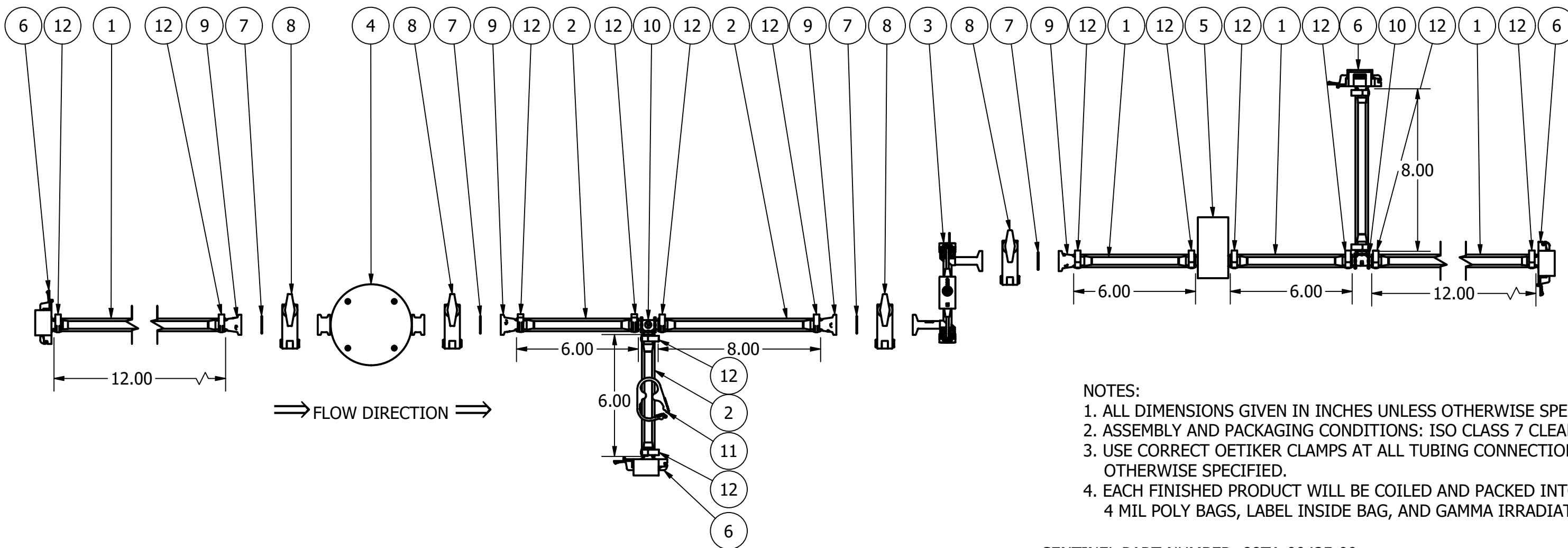

### NOTES:

1. ALL DIMENSIONS GIVEN IN INCHES UNLESS OTHERWISE SPECIFIED.
2. ASSEMBLY AND PACKAGING CONDITIONS: ISO CLASS 7 CLEANROOM.
3. USE CORRECT OETIKER CLAMPS AT ALL TUBING CONNECTIONS; UNLESS OTHERWISE SPECIFIED.
4. EACH FINISHED PRODUCT WILL BE COILED AND PACKED INTO TWO 4 MIL POLY BAGS, LABEL INSIDE BAG, AND GAMMA IRRADIATED.

SENTINEL PART NUMBER: 88TA-00425-00

REV	DATE	DESCRIPTION	DRW & CHK DATE
0			
REVISION HISTORY			

THIS DRAWING CONTAINS INFORMATION PROPRIETARY TO SENTINEL PROCESS SYSTEMS, INC. ANY REPRODUCTION IS EXPRESSLY PROHIBITED EXCEPT AS SENTINEL PROCESS SYSTEMS, INC. MAY AUTHORIZE IN WRITING.

UNLESS OTHERWISE SPECIFIED DIMENSIONAL TOLERANCES ARE:  
 0 - 12" ±1"  
 12 - 36" ±2"  
 36 - 72" ±3"  
 72" - UP ±5%

THIRD ANGLE PROJECTION



### RCTFF FEED LINE TUBING ASSEMBLY

DRAWN BY: TWG	DATE: 1/10/2023	SCALE: 1/5	PAGE NUMBER: 1 OF 1
DRAWING NO: <b>88TA-00425-00</b>		REVISION: <b>0</b>	



3265 Sunset Lane, Hatboro, PA 19040  
 P: 800.345.FLOW F: 888.329.9669 W: sentinelprocess.com